


Franklin County Public Safety
Application of Tourniquets (MAT™) A-6

| B | E | I | P | |
|---|---|---|---|--|
| | | | | Tourniquets are typically used as a last resort to control bleeding on extremity injuries or when severe, life threatening extremity hemorrhage exists. Commercial tourniquets, such as the Mechanical Advantage Tourniquet (MAT™), are preferred to as improvised tourniquets. The following steps should be performed in controlling extremity hemorrhage. The duration of these steps should last no longer than five (5) minutes before applying the MAT™. However, if severe, life threatening extremity exists, Intermediate and Paramedic providers may move to the use of the MAT™ immediately, to control the bleeding. |
| S | S | S | S | Direct Pressure |
| S | S | S | S | Elevation |
| S | S | S | S | Pressure Point |
| O | O | S | S | Tourniquet (MAT™) |

Patient Bleeding

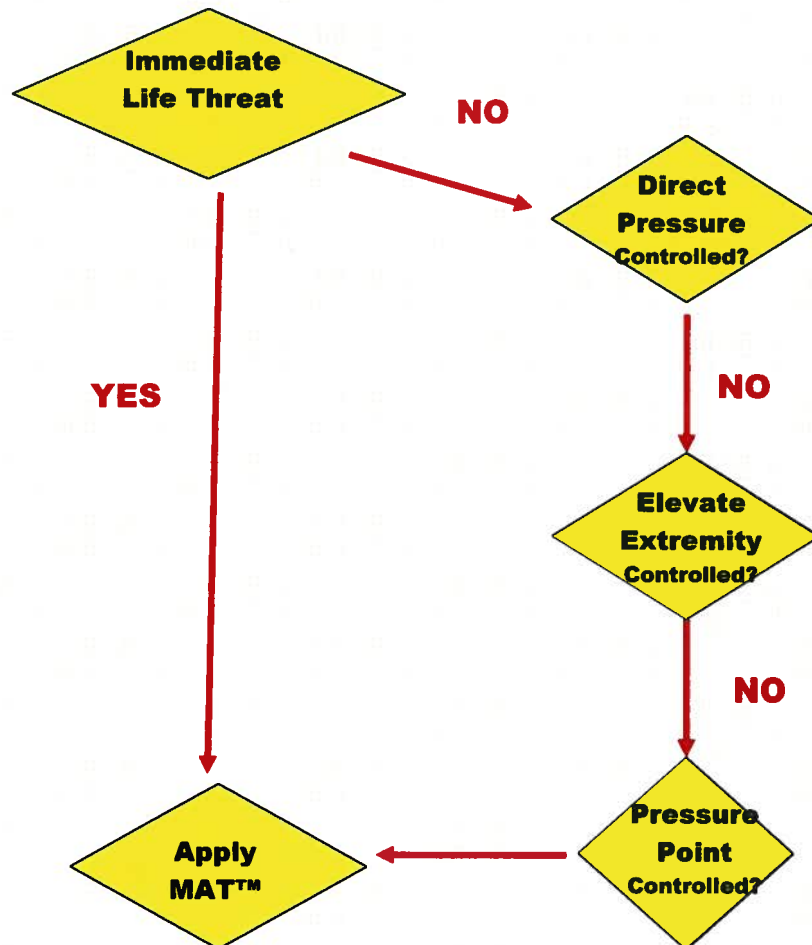


Franklin County Public Safety Application of Tourniquets (MAT™) A-6

INSTRUCTIONS AND USE

MAT[®]

Tourniquet

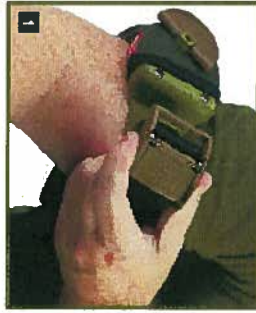

1 Place around limb, above injury.




2 "Click" buckle into place.



3 Pull strap down snugly.




4 Rotate turnkey to tighten and securely occlude blood flow.

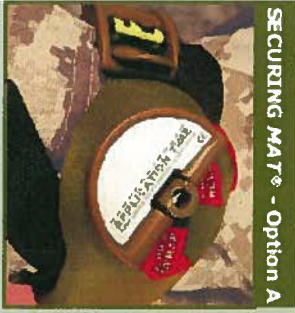


SECURING MAT[®] - Option A

Wind strap around turnkey stem and fold turnkey down.

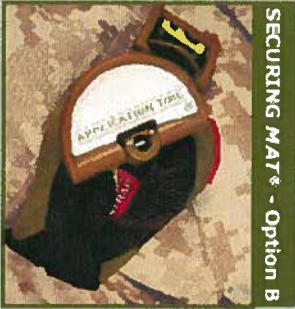

SECURING MAT[®] - Option B

Tuck strap underneath and tighten.

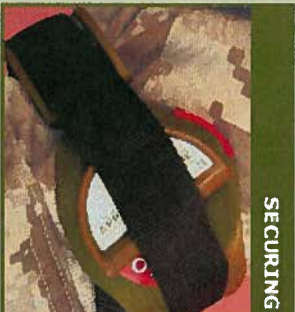


SECURING MAT[®] - Option C

Wrap excess strap over turnkey and tuck under.

Press the Release tab or lift on the buckle to release.

PYNG MEDICAL
2-12811 Cleveland Road
Reston, VA 20190
+1-800-348-9744
info@pyng.com www.pyng.com

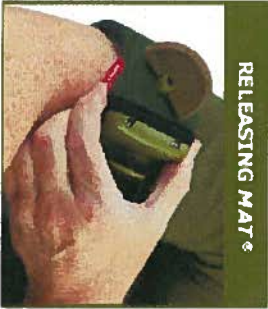

RELEASING MAT[®]

IMPORTANT
Make sure that MAT[®] is **FULLY RESET** before each application.
TO RESET
While pushing in **RED RELEASE** button, pull strap all the way out.

PM-0318